

Chemistry workshops and investigations

Make a DVD
smartphone
spectrometer

Technician Guide

To make one of these smartphone spectroscopes, each person/group will need:

A Blank DVD-R (or part of)	Scissors
Template printed on black card	Sticky tape or glue
Scalpel or craft knife	Methylated spirits and cotton wool (or a clean, soft cloth)
Black electrical tape	Phone or tablet


To use the template:

Cut along the solid lines. Use a scalpel/craft knife to cut as narrow a slit as you can in the front flap.


Fold along the dotted./dashed lines (up or down according to the instructions).

Tape your diffraction grating over the square cut-out in the 'door' – on the inside of it.

Except for the diffraction grating door, glue or tape all flaps down onto the outside.

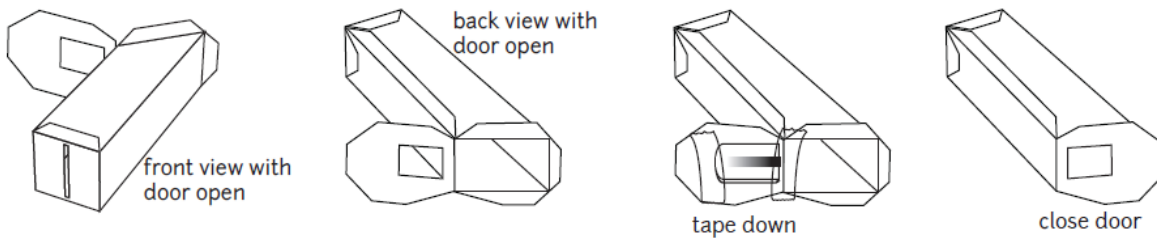

Template – print out on black card

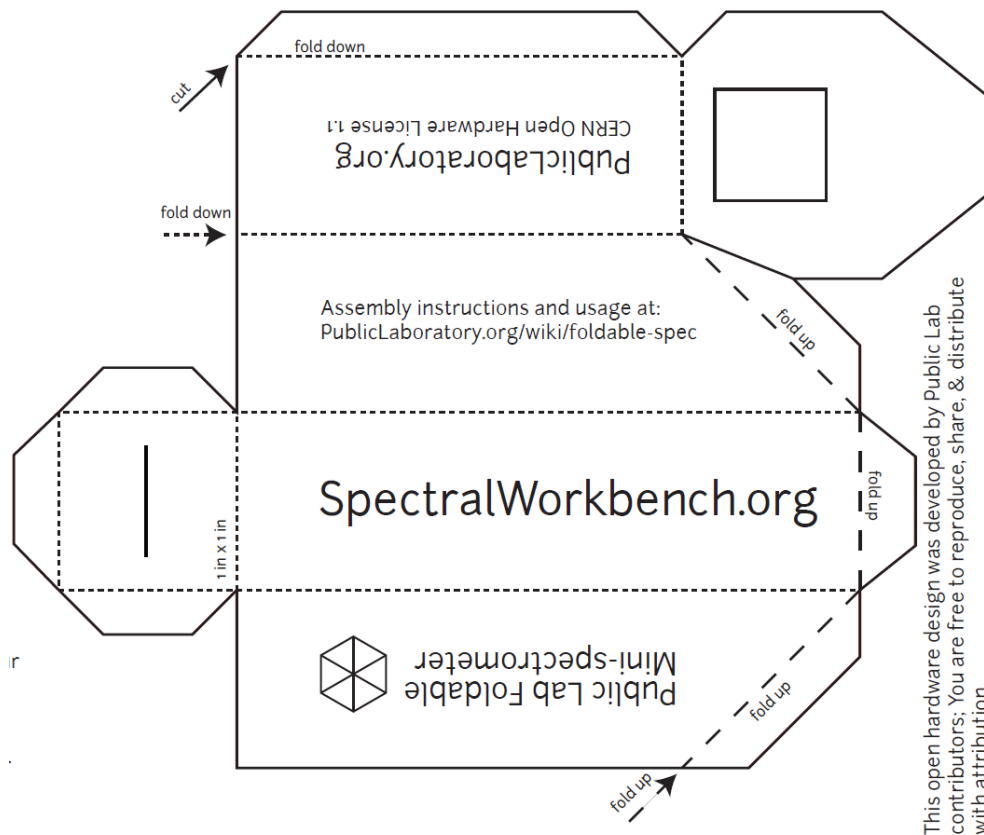

1) Preparing the diffraction grating.

You will need a recordable DVD (DVD-R).

Cut the DVD-R into halves or quarters with a pair of scissors and then you should be able to get a fingernail in and separate it into two layers: the metallic layer with the label on and a clear plastic layer with a purple sheen.

(Figure 1)

You need the clear layer. Try to handle the surface as little as possible to avoid marking it.



Figure 1: DVD layers

You can remove the purple colouring using some ethanol (methylated spirit works well) and cotton wool. (Figure 2)

Cut a piece from the edge of the now clear plastic, approximately 1cm x 1cm. This is your diffraction grating.

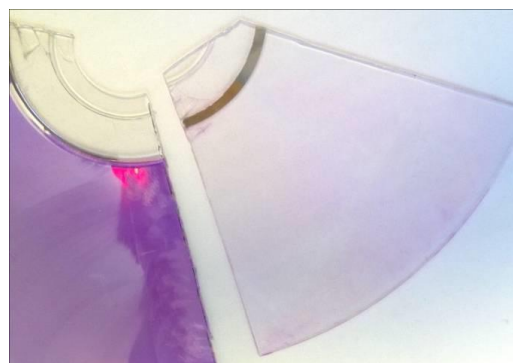


Figure 2: Cleaning the dye off

2) Getting and preparing the template

You can find the template on the publiclab website (but a copy is on the back page of this document) It refers to printing on 11" x 8" paper but A4 seems to work fine.

(<http://publiclab.org/sites/default/files/8.5x11mini-spec3.8.pdf>)

You need to print this out onto black card: card so it is stiff enough and black to reduce internal reflection. You might think printing black lines onto black card is a daft thing to do but it is still possible to see the design clearly enough to cut it out.

Cut round the outside, scissors are probably easiest.

Use a scalpel or sharp craft knife to cut out the square (top right of the diagram) where the grating goes. Then use the scalpel/knife to cut a narrow slit (lower left of the diagram) - the narrower the better.

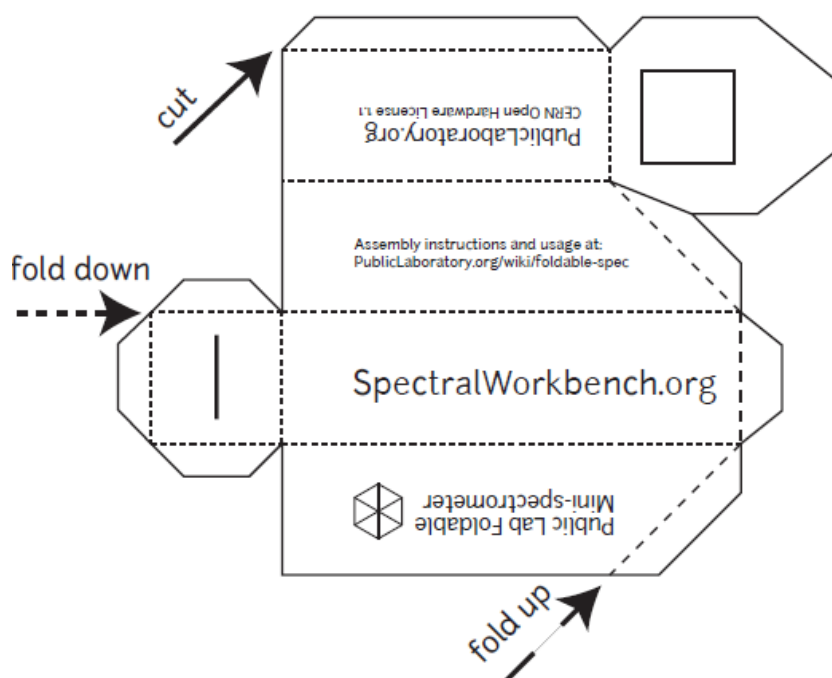


Figure 3: The template (not actual size)

3) Assembling the spectrometer

Use tape or glue to fix your 1 cm x 1 cm piece of diffraction grating on the inside of the 'door', over the square hole you cut out for it. The lines need to be horizontal so make sure that the edge of the DVD-R is at the top or bottom rather than to one side. Then fold the device as shown in the diagram and tape or glue it shut (tape is probably easier and black electrical tape is best).

You will now have a rectangular box, approximately 10 cm long with a slit at one end (Figure 4) and a slope at the other end with a window containing the diffraction grating (Figure 5)



Figure 4: The slit

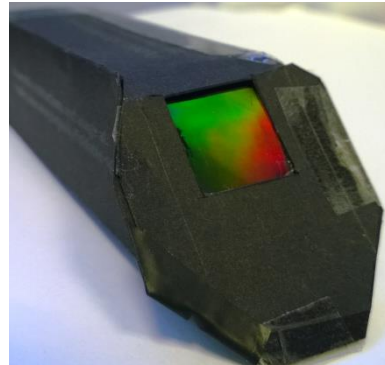


Figure 5: The diffraction grating

Cuvette adaptor

To make a cuvette adaptor, each person/group will need:

A spectroscope (made as above)	Scissors
Template printed on black card	Sticky tape or glue
Scalpel or craft knife	Black electrical tape
A cuvette	

To fit a cuvette to the front of the spectroscope.

Print the template (Figure 11) on black card.

Cut out and fold.

Solid lines are cuts, dashed lines are folds (all folds are upwards)

Use black electrical tape to hold the structure together.

Fit the adaptor over the front of a spectroscope (as shown in Figure 12). It is a lot easier if you place an empty cuvette in to make sure the fit is correct.

Secure it in place with more electrical tape

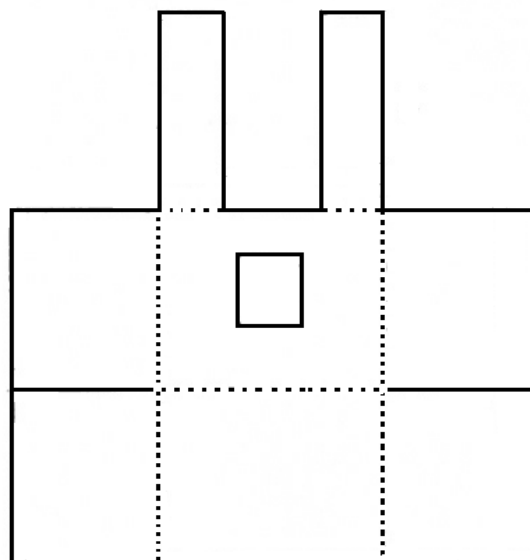


Figure 11: Cuvette adaptor template



Figure 12: Cuvette adaptor fitted



SSERC Risk Assessment (revised version November 2009)

(based on HSE '5 steps to risk assessment')

2 Pitreavie Court, South Pitreavie Business Park, Dunfermline KY11 8UB

tel : 01383 626070 fax : 01383 842793

e-mail : sts@sserc.org.uk web : www.sserc.org.uk

Activity assessed	Making a DVD Spectroscope
<i>Date of assessment</i>	8/4/2015
<i>Date of review (Step 5)</i>	
<i>School</i>	
<i>Department</i>	

Step 1	Step 2	Step 3		Step 4		
<i>List Significant hazards here:</i>	<i>Who might be harmed and how?</i>	<i>What are you already doing?</i>	<i>What further action is needed?</i>	<i>Action by whom?</i>	<i>Action by when?</i>	<i>Done</i>
Being cut during construction	Pupil using shark knife/scalpel	Follow normal procedures for using sharp implements.				
Possible damage to eyes if looking directly at too bright a light source	Pupil using spectrometer	Unlikely as spectrometer is attached to phone rather than being looked through. Warn pupils not to point it at or look directly at the sun or any other bright source				

Description of activity:

A template is cut from black card and folded into shape. A piece from a DVD-r is cut out and used as a diffraction grating. The assembled device is fixed over the lens of a mobile phone camera.

Additional comments:

Attribution

This work is based on a work by Public Lab (www.publiclab.org) The original work and this derivation is licenced under a Creative Commons Attribution-ShareAlike 3.0 Unported licence, CC BY-SA 3.0. (<http://creativecommons.org/licenses/by-sa/3.0/>)